



# SOMA MEDICAL SDN BHD


ULTRAVIOLET GERMICIDAL IRRADIATION

UVC CLEANING SYSTEMS

# UVC DEVICE M20




- **Multiple product offerings**
  - Whole Room Mobile Units – casters and handcart
- **Technology**
  - Proven UVC technology, highly reliable and powerful – 2340 watts
  - High power for effective germicidal treatment - Product efficacy - calculated and automated UVC dosage - 18 emitters
  - Pathogen inactivation - greater than a log 4 reduction of sporous pathogens on high touch surfaces
- **Product intelligence and ease of operations**
  - Proprietary based RFID tracking and web based reporting infrastructure
  - Simple operations with easy to use remote, no complicated tablets or computers
  - Safety – emitters – passive infrared AND ultrasonic
  - Efficient to operate, all from a 3 button simple remote
- **ROI and Low Cost of Ownership**
  - The average room disinfection time is less than ten (ten) minutes with linked devices (up to 8 linked devices possible).
  - Labor time savings – no need to move units around the room for multiple treatments
  - UVC will sanitize ~ 400 % more rooms in the same time over HPV products reducing the spread of infection in more locations and more often.
  - Not just focused on product but focused on environmental services and infection control process and systems SOLUTIONS.
- **Industry Validation**
  - Independent data shows UVC Cleaning product rate high in customer satisfaction and ease of use - ECRI Institute

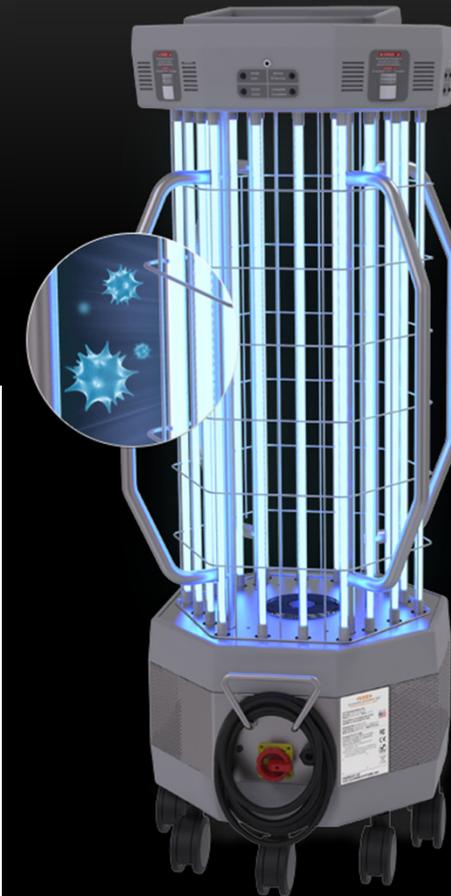

# UVC DEVICE M20 – SENSOR TECHNOLOGY

UVC Cleaning Systems uses four (4) external UVC sensors mounted on the outer surface of the device.

- These sensors measure the reflective and direct energy produced by the mobile UVC device.
- A reading is taken every one (1) second; calculates an accumulated energy threshold.
- device turns off and signals the operator “Treatment Complete” at the predetermined energy threshold level.

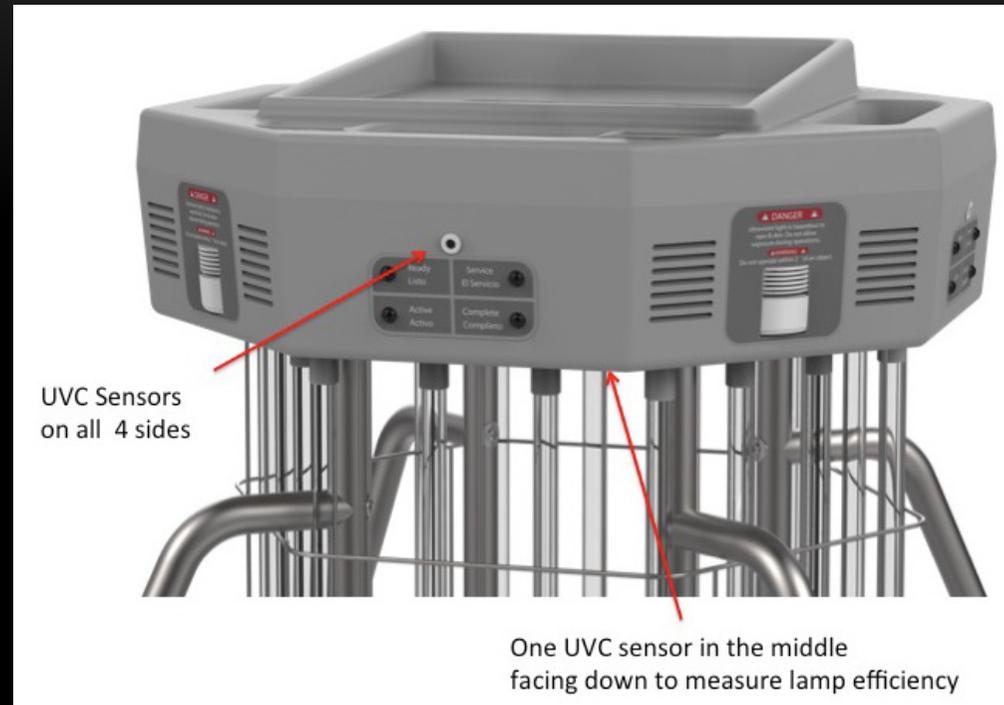


A fifth UVC sensor is mounted in center of the top assembly facing down

- measuring the amount of UVC energy produced by the number of lamps in the product.
- This sensor takes a reading every one-second measuring the efficiency of the UVC lamps.

The system is calibrated to achieve a log 4 or greater reduction of Clostridium Difficile (C. diff) on all high touch surfaces and to treat an area 23' x 23' = 529 square feet (6.7m x 6.7m = 44.89 square meters) with a single placement.

# PRODUCT APPLICATIONS

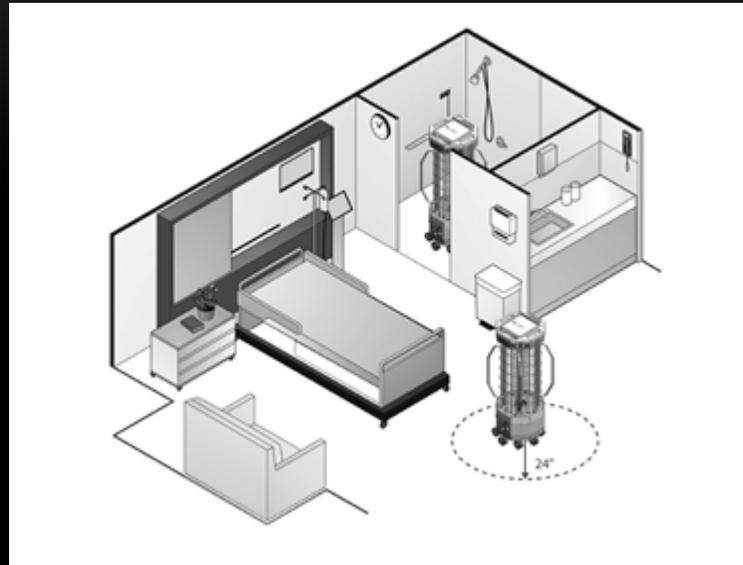

UVC Cleaning Systems has taken a different approach to create more robust designs that are effective for hospitals as well as other environments. These devices can be used in laboratories, food processing plants, hotels, cruise ships, athletic facilities, out patient surgical centers and many other environments where harmful pathogens exist. Linked or single devices provide whole room treatment where pathogens exist.

# APPLICATIONS



# TRANSPORTATION



This device is designed to easily move from facility to facility, through parking lots, and up and down stairs. This is for customers that need to be in multiple locations. One person can transport this model in a van or truck.

# OPERATION



**All UVC Cleaning Systems equipment is marked with clear, easy-to-read warning labels, and comes with a quick start guide and a comprehensive user manual to ensure safe operations.**

Machine status is indicated on all sides of the device coupled with a simple three button remote control. These simple operations minimizes confusion and ensures a safe, fast, and effective UVC treatment.

# EFFICACY



Research has proven that higher UVC energy levels reduce the population of microorganisms more than lower UVC energy in the same period of time.

# TEST REPORT

Test Facility  
Antimicrobial Test Laboratories  
1304 W. Industrial Blvd  
Round Rock, TX 78681  
(512) 310-8378

The test microorganism(s) selected for this test:

- *Staphylococcus aureus*
- *Listeria monocytogenes*
- *Clostridium difficile*

## Summary of the Procedure

1. Overnight cultures were centrifuged at 1,000 RPM (as applicable) for 10 minutes, resuspended in sterile R/O water, and supplemented with 5% FBS.
2. Stainless steel carriers (1"x3") were inoculated with 0.01 ml of the culture.
3. Inoculum was spread over approximately 10 cm<sup>2</sup> of each carrier.
4. Carriers were left to dry at room temperature for 10 minutes in Petri dishes with lids ajar.
5. Visibly dry carriers/test microorganism were harvested in 20 ml D/E to determine the initial numbers control.
6. Test carriers were treated at a specified distance and time period then harvested in 20 ml D/E.
7. Standard dilution and pour plating techniques were used for all enumerations.
8. Calculations are based off of the averaged initial and final control numbers compared to CFU/carrier recovered from treated test carriers.

# Product Background – what we inactivate

## Why we are the BEST choice

### Bacteria:

Bacillus anthracis - Anthrax  
Bacillus anthracis spores - Anthrax spores  
Bacillus magaterium sp. (spores)  
Bacillus magaterium sp. (veg.)  
Bacillus paratyphus  
Bacillus subtilis spores  
Bacillus subtilis  
Clostridium Difficile  
Clostridium tetani  
Corynebacterium diphtheriae  
Ebertelia typhosa  
Escherichia coli  
Listeria  
Leptospira canicola - infectious Jaundice  
Micrococcus candidus  
Micrococcus sphaeroides  
Mycobacterium tuberculosis  
Neisseria catarrhalis  
Phytomonas tumefaciens  
Proteus vulgaris

### Bacteria:

Pseudomonas aeruginosa  
Pseudomonas fluorescens  
Salmonella enteritidis  
Salmonella paratyphi - Enteric fever  
Salmonella typhosa - Typhoid fever  
Salmonella typhimurium  
Sarcina lutea  
Serratia marcescens  
Shigella dysenteriae - Dysentery  
Shigella flexneri - Dysentery  
Shigella paradysenteriae  
Spirillum rubrum  
Staphylococcus albus  
Staphylococcus aureus  
Staphylococcus hemolyticus  
Staphylococcus lactis  
Streptococcus viridans  
Vibrio comma - Cholera

### Viruses:





Bacteriophage - E. Coli  
Infectious Hepatitis  
Influenza  
Poliovirus - Poliomyelitis  
Tobacco mosaic  
Norovirus

### Yeast

Brewers yeast  
Common yeast cake  
Saccharomyces cerevisiae  
Saccharomyces ellipsoideus  
Saccharomyces spores

### Molds:

Aspergillus flavus  
Aspergillus glaucus  
Aspergillus niger  
Mucor racemosus A  
Mucor racemosus B  
Oospora lactis  
Penicillium expansum  
Penicillium roqueforti  
Penicillium digitatum  
Rhizopus nigricans

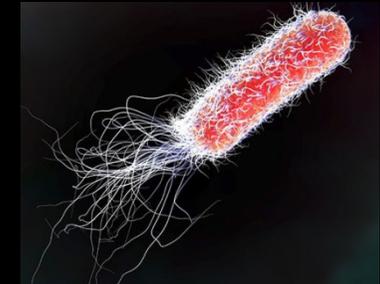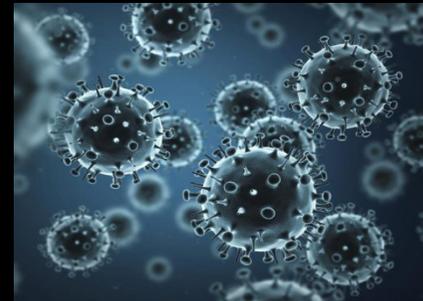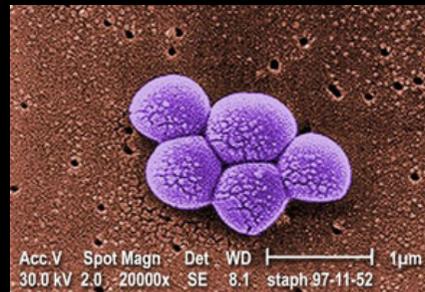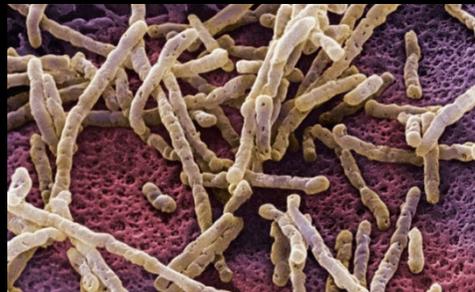

# TEST REPORT – BACTERIA TESTED



## ***Staphylococcus aureus* 6538**

This bacterium is a Gram-positive, spherical-shaped, facultative anaerobe. Staphylococcus species are known to demonstrate resistance to antibiotics such as methicillin. *S. aureus* pathogenicity can range from commensal skin colonization to more severe diseases such as pneumonia and toxic shock syndrome (TSS). *S. aureus* is commonly used in several test methods as a model for gram positive bacteria. It can be difficult to disinfect but does demonstrate susceptibility to low level disinfectants.



## ***Listeria monocytogenes***

This bacteria is a Gram-positive, rod shaped, facultative anaerobe that is motile due to the presence of flagella. These bacteria are common cause of the foodborne illness listeriosis, which can be fatal. Listeriosis can cause meningitis and sepsis and is particularly dangerous to pregnant women and unborn infants. *Listeria monocytogenes* is pervasive and can be found in soil, water, and certain livestock animals. They can resist both heat and freezing and can survive for several years.



## ***Clostridium difficile* 43598**

This bacteria is a Gram-positive, rod shaped, endospore generating obligate anaerobe. Clostridium species are part of the normal human gut flora that produce spores which are highly resistant to chemical and environmental conditions. *C. diff* is commonly associated with hospital acquired infections and is known to cause antibiotic assisted colitis. Because of its high resistance to antimicrobials, *C. difficile* is a benchmark bacteria for sporicidal and sterilant activity of chemicals.


# TEST REPORT - RESULT

## *Staphylococcus aureus*

Species	Contact Time	Contact distance	% reduction
<i>S. Aureus</i> (run 1)	10 minutes	3 meters	99.9987%
		5 meters	99.9967%
		5 meters perpendicular	99.9975%
<i>S. Aureus</i> (run 2)	18 minutes	3 meters	99.9998%
		5 meters	99.9989%
		5 meters perpendicular	99.9993%

# TEST REPORT - RESULT


## *Staphylococcus aureus*

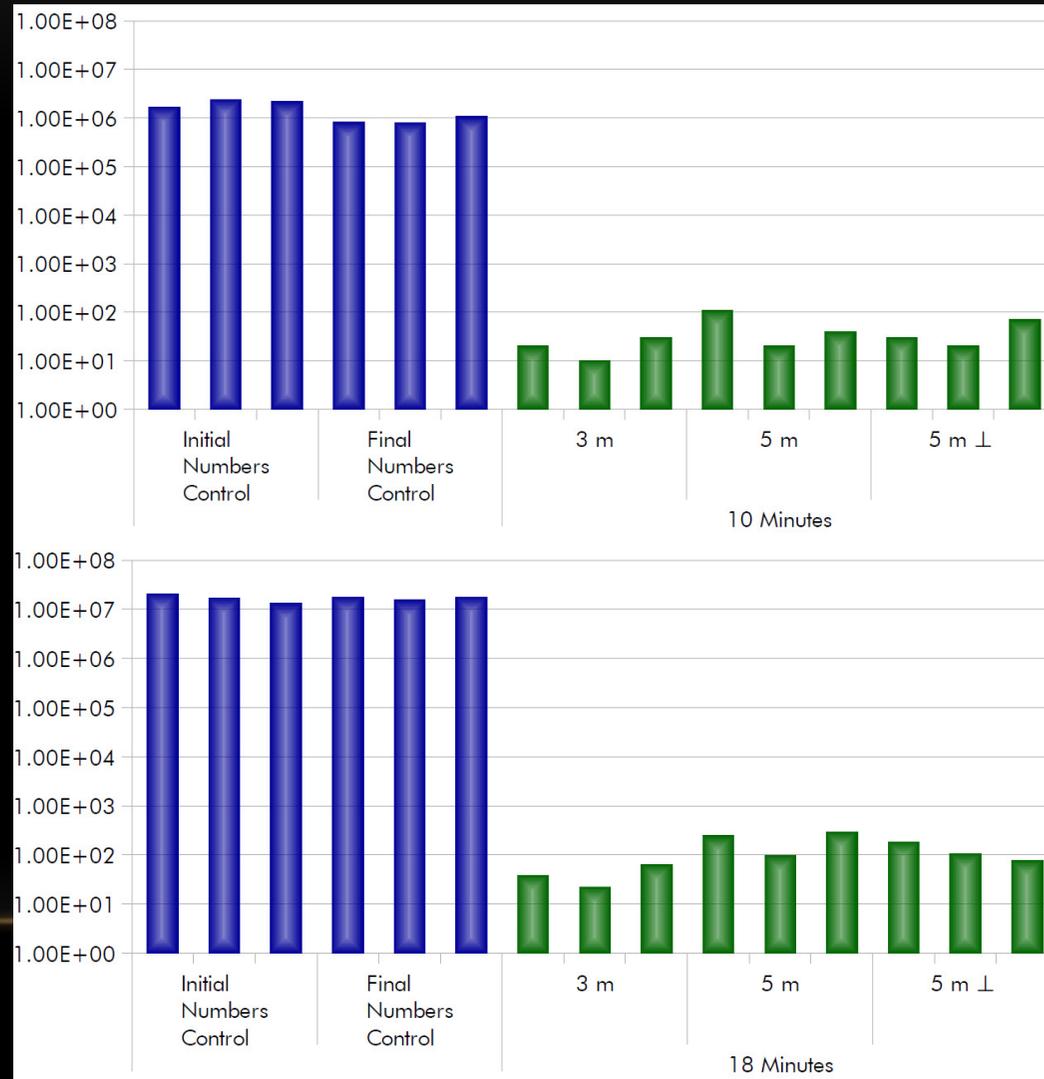

# TEST REPORT - RESULT

## *Listeria monocytogenes*



Species	Contact Time	Contact distance	% reduction
<i>L. monocytogenes</i> (run 1)	18 minutes	5 meters	99.9998%
		5 meters perpendicular	99.978%
<i>L. monocytogenes</i> (run 2)	25 minutes	5 meters	99.9988%
		5 meters perpendicular	>99.9997%

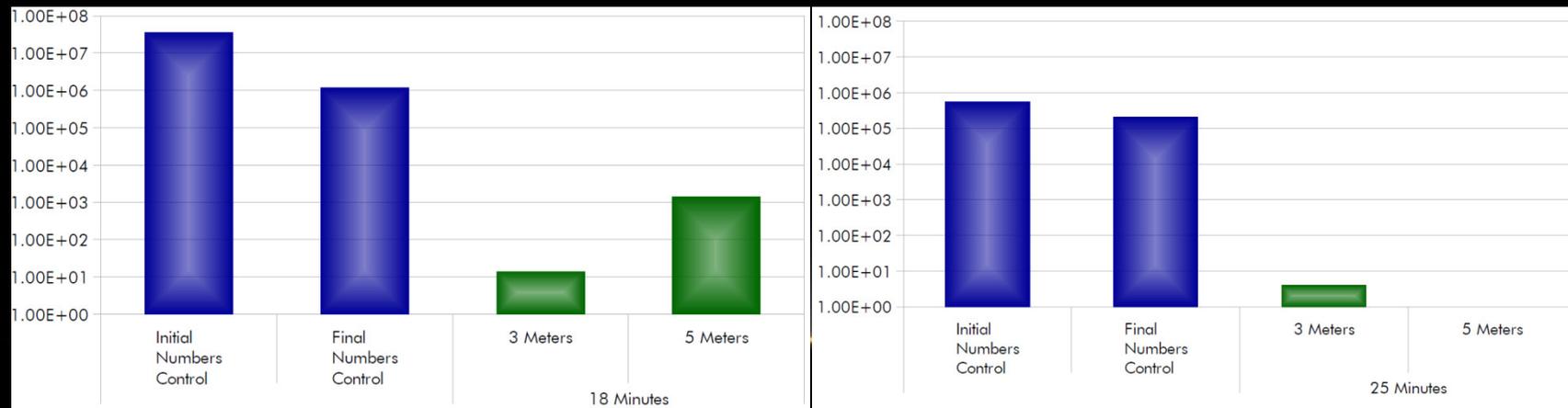


# TEST REPORT - RESULT

## *Clostridium difficile*

Species	Contact Time	Contact distance	% reduction
<i>C. difficile</i> (run 1)	12 minutes	3 meters	99.03%
		5 meters	99.16%
<i>C. difficile</i> (run 2)	55 minutes	3 meters	99.994%
		5 meters	99.992%

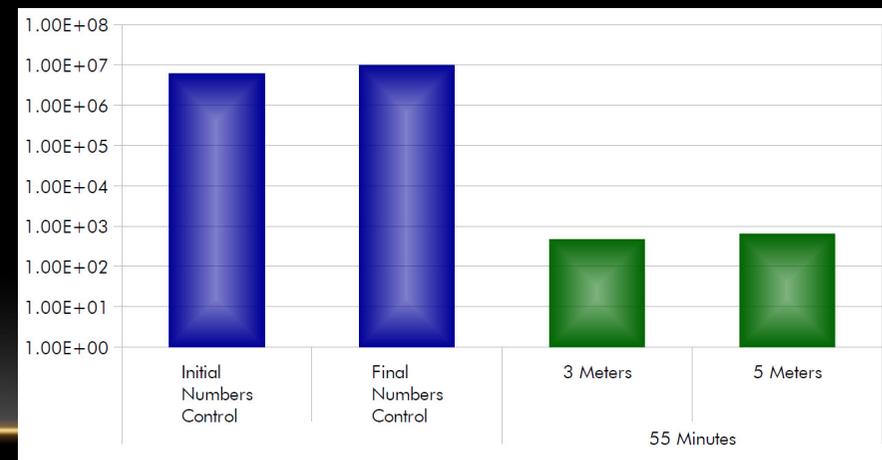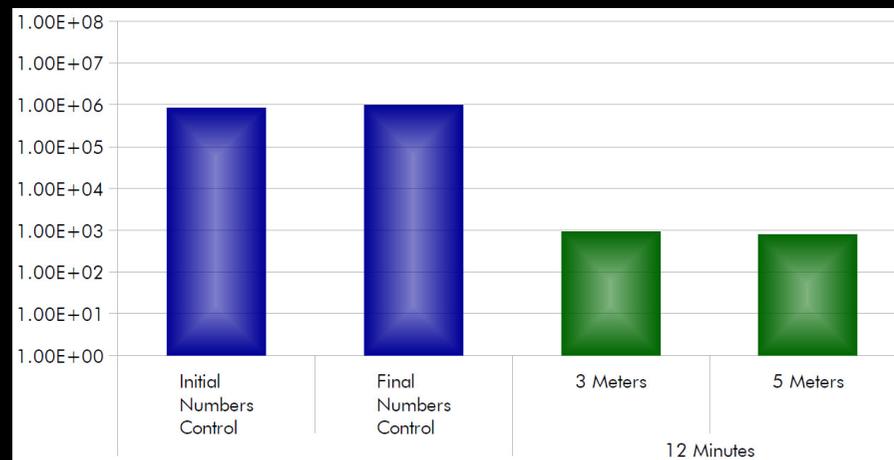

# SAFETY

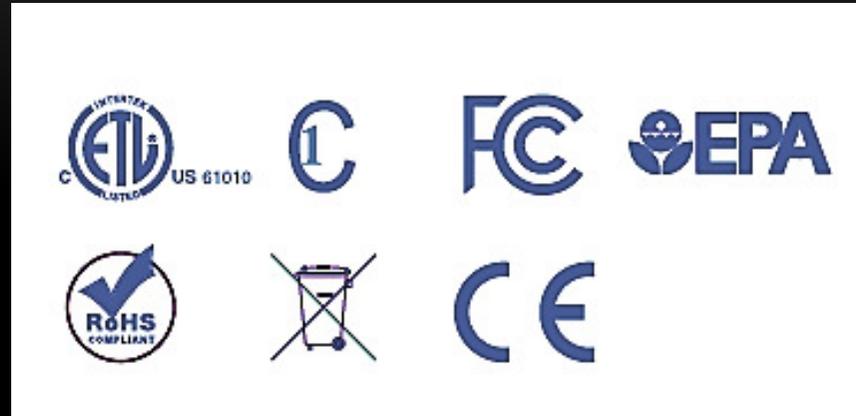

UVC Cleaning Systems is committed to producing the safest, most effective UVC disinfection devices on the market. The equipment is ETL, CE, CSA, and FCC certified and is WEEE, RoHS, and EPA compliant.

# SM NANO 1152 TITANIUM DIOXIDE



[Physicochemical Data Sheet]	
Product Series	Nano TiO <sub>2</sub> Sol Coating Agent (SM1152)
Appearance	Transparent liquid
Dispersive type	Solution
Odor	None
PH	7-8.5
Boiling Point	100°C/212°F
Volatility	None
Freezing Point	0°C/32°F
Flash Point	Non flammable
Average primary particle size - Acc. to GB/T 19591-2004	< 4nm
Crystal structure - Acc. to GB/T 19591-2004	Anatase
Specific surface area (BET) - Acc. to ISO 9277:1995	160± 30m <sup>2</sup> /g
Coagulation Index - Acc. to GB/T 19591-2004	2-4
Material academic duration	Permanent
Coating duration - Acc. to outdoor simulation environment	> 2 years
Primary drying time	30 minutes
Final setting time	2 weeks
Saturated stream pressure	2333Pa acc. to H <sub>2</sub> O 1 PN 20°C
Opposite stream density	< 1.0 acc. to H <sub>2</sub> O
Solubility	Dissolve in water, miscible in oil
True specific gravity	1.0075 – 1.01
Viscosity, dynamic	1.0050 mPa.s
Vaporize velocity	< 1.00 acc. to H <sub>2</sub> O

# ULV INTELLIGENT ELECTRIC COLD FOGGER

## Features

- Good atomization effect with uniformed droplet.
- With its high efficiency and pressure atomization control, the sprayer saves on chemical and any kind of leakage is prevented.
- Easy to carry and operate.
- Stable performance.

## Parameter

Motor	220V/50Hz
Power	800W
Tank capacity	2 Litres
Effective range	6 Metres
Noise factor	70dB
Particle size	0.4 to 1.8 Microns
Net weight	3.1 kg
Dimension	581 x 150 x 229.5 mm



# ALS TEST REPORT ON TITANIUM DIOXIDE

## ALS TECHNICHEM (M) SDN BHD

(117964-P)

9, Jalan Astaka U8/84, Seksyen U8, Bukit Jelutong, 40150 Shah Alam, Selangor.

Tel: (603) 7845 8257 Fax: (603) 7845 8258 E-mail: info@alsmalaysia.com



ALS Technichem

### CERTIFICATE OF ANALYSIS

#### MICROBIAL CHALLENGE TEST

Day Organism	Inoculum Count (cfu/ml)	Count Recovered (cfu/ml)	% Reduction
<i>Staphylococcus aureus</i> ATCC 6583	$4.5 \times 10^6$	$2.0 \times 10^4$	99.56
<i>Pseudomonas aeruginosa</i> ATCC 9027	$6.2 \times 10^6$	$3.0 \times 10^4$	99.52
<i>Escherichia coli</i> ATCC 8739	$7.3 \times 10^6$	$3.7 \times 10^4$	99.49
<i>Candida albicans</i> ATCC 10231	$4.5 \times 10^5$	$3.9 \times 10^3$	99.13
<i>Aspergillus niger</i> ATCC 16404	$3.7 \times 10^5$	$3.3 \times 10^3$	99.11
<b>Method References</b>	<b>USP &lt;51&gt;</b>		

Remark: 1.0 ml inoculum from  $10^6$  was injected in 100ml sample

Contact time: 1 Hour

Calculation for % reduction:  $\frac{\text{Inoculum Count} - \text{Count Recovered}}{\text{Inoculum Count}} \times 100$

USP: United State Pharmacopoeia

# THANK YOU



*102, Lorong Maarof,  
Bangsar, 59000 Kuala Lumpur  
Malaysia*

*T: +603 2287 5790 F: +603 2287 6790*

*[www.somamedical.net](http://www.somamedical.net)*

*Email: [info@somamedical.net](mailto:info@somamedical.net)*